

Tutorial 3 – Muscle-tendon Personalization

The Muscle–tendon Model Personalization tool finds an optimal set of subject-specific muscle–tendon properties and muscle activations from EMG, joint kinematic, and joint moment data by balancing optimization cost function terms related to muscle properties, similarity of properties among grouped muscles, and matching of EMG-driven and experimental inverse dynamics joint moments. Muscle activation and force predictions are sensitive to optimal muscle fiber length and tendon slack length. Therefore, reliable personalization of these parameters is essential for generating reliable predictions of muscle activations and forces during predictive simulations of movement.


The inputs to the MTP tool are a post-JMP OpenSim model as well as IK motion, ID load, muscle–tendon length and velocity, and muscle moment arm data from one or more motion trials of interest.

Before running MTP:

1. Open the OpenSim model **UF_Subject_3_reduced_muscles.osim** in the OpenSim GUI.
2. Under the *Forces* tab on the model, explore the muscles available.
3. Take note of the extra groups added.
 - a. These are added for organization so that MTP/NCP knows which model muscles to group together in the optimization.
 - b. The four important groups are:
 - i. Activation Muscle Groups – Muscles that we would expect to have similar activation profiles (ie lateral hamstrings; BFSH and BFLH will have similar activations to each other)
 - ii. Normalized Fiber Length Muscle Groups – Muscles that we would expect to have similar normalized fiber lengths.
 - iii. Collected EMG Muscle Groups – Muscle groups that we **do have** experimental EMG data for. These must have the same name as the respective EMG channel name your EMG data file (**preprocessed\EMGData\gait_1.sto**)
 - iv. Missing EMG Muscle Groups – Muscle groups that we **do not have** experimental EMG data for.
 - c. These groups need to be created manually in a text editor such as Notepad++, but we have example groups for lower limb models in the *NMSM Article* download on SimTK: https://simtk.org/frs/?group_id=2397

Setting up an MTP settings file:

1. Activate the NMSM GUI in OpenSim by navigating to *Tools>User Plugins*, and click **rcnlPlugin.dll**.
2. With **UF_Subject_3_reduced_muscles.osim** selected in the OpenSim GUI, navigate to *Tools>Model Personalization>Muscle-tendon Personalization*.
 - a. The following window should be opened:

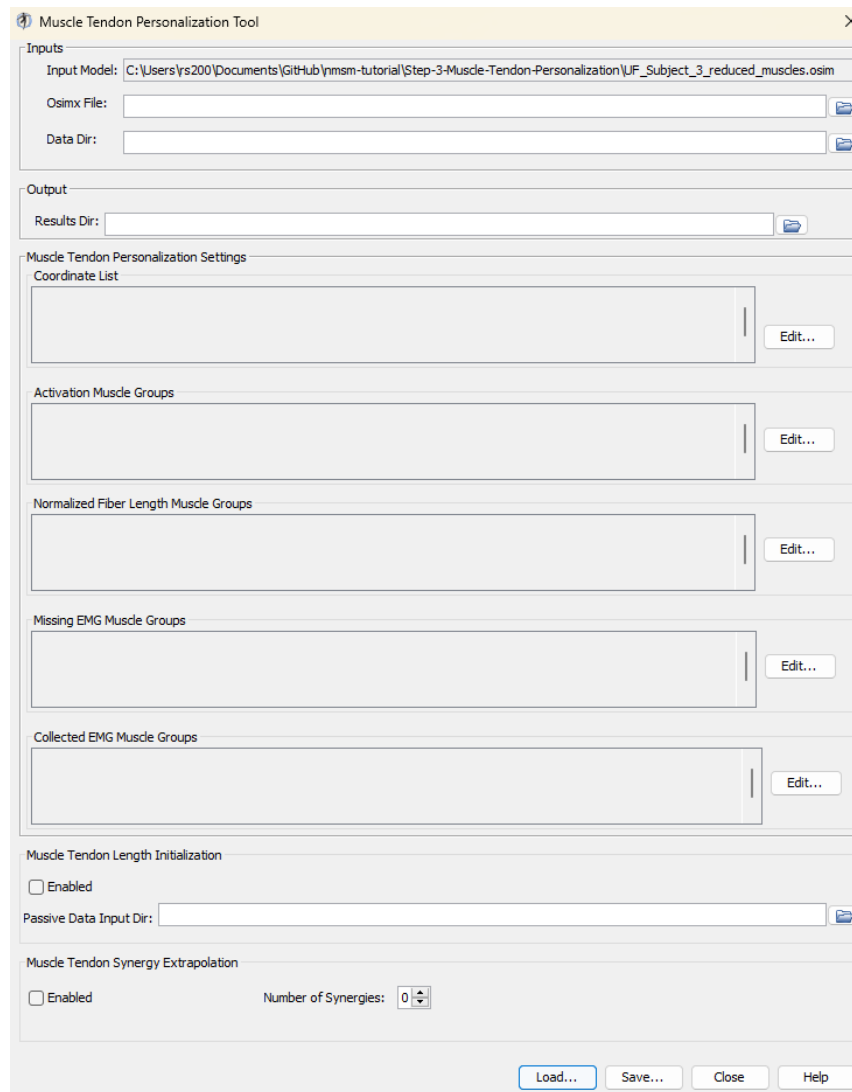

The screenshot shows the 'Muscle Tendon Personalization Tool' window. It has a title bar with a close button. The window is divided into several sections:

- Inputs:** Contains three fields: 'Input Model:' with a file path, 'Osimx File:' with an empty text box and a folder icon, and 'Data Dir:' with an empty text box and a folder icon.
- Output:** Contains one field: 'Results Dir:' with an empty text box and a folder icon.
- Muscle Tendon Personalization Settings:** This section contains five sub-sections, each with a list box and an 'Edit...' button:
 - Coordinate List:** An empty list box.
 - Activation Muscle Groups:** An empty list box.
 - Normalized Fiber Length Muscle Groups:** An empty list box.
 - Missing EMG Muscle Groups:** An empty list box.
 - Collected EMG Muscle Groups:** An empty list box.
- Muscle Tendon Length Initialization:** Contains a checkbox labeled 'Enabled' (which is unchecked) and a 'Passive Data Input Dir:' field with an empty text box and a folder icon.
- Muscle Tendon Synergy Extrapolation:** Contains a checkbox labeled 'Enabled' (which is unchecked) and a 'Number of Synergies:' field with a spinner box set to 0.

At the bottom of the window are four buttons: 'Load...', 'Save...', 'Close', and 'Help'.

3. Leave the *Osimx File* field empty. This tool outputs an Osimx file, but we do not have one to work with yet. If this field is filled out, the MTP tool will concatenate new elements to the existing Osimx file.
4. Set the *data directory* to be **preprocessed**.
5. Set the *results directory* to be **MTPResultsV1**.
6. For the *coordinate list*, select: (**hip_flexion_r**, **knee_angle_r**, **ankle_angle_r**).

7. For *activation muscle groups*, select (**HipFlexorsActivationGroupR**, **GlutmaxActivationGroupR**, **HamslatActivationGroupR**, **VasActivationGroupR**, **GasActivationGroupR**).
 - a. Tip: The filter box at the top is very helpful to filter only activation groups.
8. For *normalized fiber length groups*, select (**GlutmaxNormalizedFiberLengthGroupR**, **HamsNormalizedFiberLengthGroupR**, **VasNormalizedFiberLengthGroupR**, **GasNormalizedFiberLengthGroupR**).
9. For *missing EMG muscle groups*, select (**HipFlexorsMissingEMGChannelGroup**)
10. For *collected EMG muscle groups*, select (**GlutMaxLat**, **RecFem**, **BicFemLong**, **BicFemShort**, **VastMed**, **GasMed**, **Sol**, **TibAnt**).
11. **Enable Muscle Tendon Length Initialization** (MTLI).
12. Set the *passive data input directory* to **passive_moment_data**.
13. **Enable Muscle tendon Synergy Extrapolation** (SynX) with **3 synergies**.
14. Save this settings file as **MTPSettingsV1.xml**.
15. Open **MTPSettingsV1.xml** in a text editor of your choice and explore the settings file.

Running MTP:

1. Open MATLAB and open **runMTP.m** in your tutorial directory.
2. Open the project file (**Project.prj** inside your installation of nmsm-core).
3. Run the MATLAB section labelled **Run MTP V1**.
 - a. With the section selected, press shift+enter to run a section.

Post MTP Analysis:

1. Look through the plots created by the script. If everything was done correctly, there should be 6 plots.
 - a. Plot 1 – Joint Moment Matching: Joint moments generated by muscle forces (with and without SynX) compared to Inverse Dynamics joint moments.
 - b. Plot 2 – Muscle Activations: Muscle excitations and activations (with and without SynX) for all muscles included in the MTP run.
 - c. Plot 3 – Normalized Fiber Lengths: Normalized fiber lengths for all muscles included in the MTP run. Red dashed lines at 1.0 and 0.6 indicate the optimal working range for muscles.
 - d. Plot 4 – Passive Muscle Force: The passive force generated by each muscle in the MTP run. This plot is only generated if MTLI is enabled.
 - e. Plot 5 – Passive Joint Moment Matching: The passive moment matching achieved by the MTLI optimization.

- f. Plot 6 – Hill-type Muscle Model Parameters: The muscle model parameters for all optimized muscles in this MTP run. These parameters are the primary output of the MTP tool.
2. Using a text editor, explore **UF_Subject_3_reduced_muscles_mtp.osimx** in **MTPResultsV1**. This file is used as an input for the Neural Control Personalization (NCP) tool, and for Treatment Optimization tools if using synergy controls.

Experiment with different numbers of synergies:

1. Create a copy of **MTPSettingsV1.xml** and name it **MTPSettingsV2.xml**.
2. Open **MTPSettingsV2.xml** in a text editor of your choice.
3. Change `<results_directory>` to **MTPResultsV2**.
4. Inside `<MTPSynergyExtrapolation>`, change `<number_of_synergies>` to **4**.
5. Run the MATLAB section labelled **Run MTP V2**.

Run MTP without SynX:

1. Create a copy of **MTPSettingsV1.xml** and name it **MTPSettingsV3.xml**
2. Open **MTPSettingsV3.xml** in a text editor of your choice.
3. Change the `<results_directory>` to **MTPResultsV3**.
4. Inside `<MTPSynergyExtrapolation>`, change `<is_enabled>` to **false**.
5. Save this settings file as **MTPSettingsV3.xml**.
6. Run the MATLAB section labelled **Run MTP V3**.